Evolsin® Anti Stielwarze

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Rev 1 / vom 15 Januar 2025

Evolsin® Anti Skin Tags

EVOL SIN®ANTI SKIN TAGS is a medical device for self-application, designed for the simple, safe and effective removal of skin tags.

IDENTIFYING SKIN TAGS:

Skin tags are also known as soft fibromas. These are benign growths of the connective tissue. Skin tags can be identified by the following characteristics:

- Skin tags are small, stalk-shaped skin growths. They may have a thickened end resembling a small pouch.
- Skin tags are always soft. They can be touched without causing pain. Because they are soft and stalk-like, they can be gently twisted and moved with a finger.

PURPOSE: Removal of soft skin tags (fibromas).

INDICATIONS: Skin-coloured skin tags with a base width of up to 2 mm.

CONTRAINDICATIONS:

The product must **not** be used:

- For the removal of skin tags on eyelids, near the eyes, on mucous membranes or in the genital or anal areas
- To remove bleeding skin tags or skin tags showing signs of skin disease or infection
- On moles, birthmarks or skin tags with hair growth
- On irritated skin, infected or reddened areas, or skin with itching or swelling
- On individuals with skin cancer (current or past)
- On those with blood disorders or impaired wound healing
- On those taking blood-thinning medications
- On individuals with known allergies to ABS (acrylonitrile-butadiene-styrene) or silicone
- In combination with any other skin tag removal methods

PATIENT AND USER GROUP:

Adults aged 18 and over, including pregnant and breastfeeding women.

SIDE EFFECTS:

A slight tingling sensation may occur during treatment. Typically, there is little to no scarring.

SAFETY PRECAUTIONS:

Keep out of reach of children.

WARNINGS:

EVOL SIN®ANTI SKIN TAGS should only be used for the types of skin tags specified above (soft fibromas). If you are unsure whether something is, indeed, a skin tag, consult your doctor. Using the product for other purposes may negatively impact your health.

Silicone rings are intended for single use only. Do not reuse silicone rings. Dispose of them after use.

NOTE FOR PATIENTS / USER:

Any serious incident that has occurred in relation to the device must be reported to the manufacturer and the relevant authority in your EU Member State.

HOW DOES EVOL SIN®ANTI SKIN TAGS WORKS?

During use, a snug silicone ring is placed around the base of the skin tag. This cuts off its blood supply. By restricting blood flow, EVOL SIN®ANTI SKIN TAGS causes the skin tag to die off and detach, typically within a few days. New skin forms beneath the base during this process.

HOW TO USE EVOL SIN®ANTI SKIN TAGS:

Evolsin® Anti Skin Tags consists of the following components:

- 1 small application funnel for small silicone rings
- 1 large application funnel for large silicone rings
- 1 applicator rail
- 20 small silicone rings
- 20 large silicone rings
- Instructions for use

Items to prepare:

- Mild soap or disinfectant for cleansing the skin area
- Plaster to cover the treated skin area
- A mirror, if required

INSTRUCTIONS FOR USE:

For external use only. Use only as directed.

- 1) Clean your hands, the skin tag and the surrounding skin. Use a mild soap and water or a commercially available skin disinfectant. Avoid oil-based soaps or impregnated wet wipes, as residue may prevent the silicone ring from adhering properly. Dry the skin thoroughly. (Fig. 1)
- 2) Select the appropriate application funnel: Use the small funnel for small skin tags. It should fit completely over the skin tag. If the skin tag is too large for the small funnel, use the large funnel.
- 3) Choose a silicone ring. The size of the silicone ring corresponds to the funnel size. Slide the silicone ring over the tip of the application funnel to just before the wide end. (Fig. 2)
- 4) Place the application funnel over the skin tag. Ensure that the wide, circular end of the funnel is fully in contact with the skin. (Fig. 3)
- 5) Hold the applicator rail in the other hand. Slide the matching opening over the tip of the application funnel. Gently press the rail down until the silicone ring slides over the edge of the funnel and tightly encloses the base of the skin tag. (Fig. 4)
- 6) Cover the skin tag with a commercial plaster.

Fig. 1

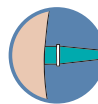
Fig. 2

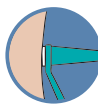
Fig. 3

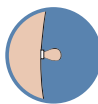
Fig. 4

IN THE FOLLOWING DAYS:

Take care during personal hygiene to avoid accidentally dislodging the silicone ring. It may take several days for the skin tag to detach.

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